SANGLES

NAUSICAA
MEDICAL

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

Verticonfort riser sling

Description

The strap is positioned at the bottom of the back (at the lumbar level) in a seated position. The strap must be taut and attached to the first or second notch (green or yellow) on the arm of the standing frame. It is essential to choose the correct size for the patient's build.

Sizes

	S	M	L	XL
A	81	94	107	123
B	21	25	29	29

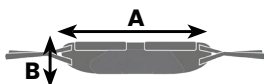SVERTICONF-PF

V03.07/26

Recommendation

The use of the device requires the presence and vigilance of a third person. This medical device must be implemented or trained by a health professional.

The strap is positioned at the bottom of the back (at the lumbar level) in a seated position. The strap must be taut and attached to the first or second notch (green or yellow) on the arm of the standing frame.

It is essential to choose the correct size for the patient's build.

Composition

- 100 % polyester jersey
- 40 mm polyester strap - 1450 DaN
- 50 kg/m³ PE foam
- 100 kg/m³ PE foam
- 20 kg/m³ polyurethane foam
- 18 kg/m³ polyurethane foam
- Coating polyurethane on polyester
- 50 mm polyester bias
- 20 mm polyester bias
- Self grip 100 % polyester

Set up

This device is intended for stand up lifts allowing to transfer a person in standing position. Before each use, it is important to check the presence of the label, the stitches of the straps, that they do not fray or break, that the device does not show any tears, snags, holes or signs of fraying (fabric and straps), that the buckles are not cracked or have no weak fasteners

Position the bottom of the backrest around the body size in order that the ears of the attachment system are on either side of the person.

Hook the two fasteners of the backrest to the stand-up lift hooks.

SETTINGS

The fastening system, with multiple hangers, allows to adjust the position of the person.

Important recommendations

- In order to get a maximum efficiency of these sort of product, it is necessary :
 - to choose the appropriate size for the patient
 - to get the best possible adjustment to the patient
- These products must not be in direct contact with a wounded skin.

Caution

Read the instruction manual before using the slings.

The backrest hangers must be hung on the same color.

Never change the attachment setting of the ears to the stand-up lift hooks, when someone is installed in the sling.

Before lifting a person with the stand-up lift, make sure the sling is securely hooked to the hooks.

End-of-life management