www.nausicaa-medical.com

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www.nausicaa-medical.com

This device is intended lifts and ceiling-mounted lifts for transferring a person to a sitting position.

It is compatible with lifts equipped with a smooth, open or closed hook system, being able to receive a 40 mm wide sling without impacting the strength of the sling on beams ranging from 2 to 4 points, which excludes hooking with a metal plate with holes.

QUICK ADJUST U SHAPE SLING

Compatible with
removable headrest
(TAMOV).

- CE certified and manufactured in France, the Fast Adjustable Comfort Sling is recommended for bed to chair transfers requiring a management of the head support.
- We recommend using it for bed to chair transfers requiring a management of the head support (with the addition of removable headrest).

Before each use, it is important to check: the presence of the label, the stitching of the straps, that they do not fray or come undone, that the device has no tears, snags, holes or signs of fraying (fabric and attachment straps), that the buckles are not cracked or have weak attachments.

Remove the aluminium bar before putting the sling in the washing machine.

SST-PF / SML-PF / SSTBA-PF / SMLBA-PF

The use of the device requires the presence and vigilance of a third person. When used by an agitated or disoriented person, monitoring becomes imperative.

This medical device must be implemented or trained by a health professional.

2424 70 V04.07/26

Important Recommendations :

- In order to get a maximum efficiency of these sort of product, it is necessary :
 - to choose the appropriate size for the patient
 - to get the best possible adjustment to the patient
- These products must not be in direct contact with a wounded skin.


QUICK ADJUST U SHAPE SLING

SIZES & DIMENSIONS (in cm)

	A	B	C
S	57	62	54
M	61	67	65
L	65	69	76
XL	70	71	86

Composition
100% Polyester Fabric
150Kg/M3 PE foam
Sling 100% Polyester

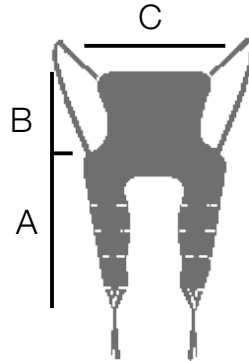CAUTION

Read the instructions for use before using the sling.
Hang to the spreader bar the loop which are at the same level.
Do not change any adjustment while the patient is hanging in the sling.
Before lifting a patient, ensure that the sling is correctly hooked on the spreader bar.

End-of-life management