NAUSICAA
MEDICAL

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Made in France
since 1992

www.nausicaa-medical.com

This medical device is intended for patient support and positioning. It is not a restraint system. Before each use, it is important to check: the presence of the label, that the device does not show any tears, snags, holes or signs of fraying and the good condition of the hook and loop.

FASTENERS

PADDED WRIST

LIGHT WRISTS SUPPORT FOR LOW-RESTLESSNESS PATIENTS

SECURITY WRIST FASTENERS

FOR A STRONGER CONTENTION FOR VERY AGITATED PATIENTS

ANKLE FASTENERS

MAINTAINS THE PATIENT'S ANKLES, WHILE ALLOWING FREEDOM OF MOVEMENT.

⚠

Restraint is exceptional and must be limited in time (requires a medical prescription). It must be accompanied by a regular assessment of the benefits/risks and increased patient monitoring. It should only be used by staff who are trained and qualified in its use. In the event of incorrect use, there is a risk of accident to the patient.

⚠

The use of the device requires the presence and vigilance of a third person. When used by an agitated or disoriented person, monitoring becomes imperative.
This medical device must be implemented or trained by a health professional.

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CE

Important Recommandations :

• In order to get a maximum efficiency of these sort of product, it is necessary :

- to choose the appropriate size for the patient
- to get the best possible adjustment to the patient

• These products must not be in direct contact with a wounded skin.

FASTENERS

Composition :
PADDED WRIST FASTENER
Combination of polyurethan foam and jersey/ cotton polyester fabric 100 % polyester fur
100 % polyester velcro strap
100 % polyester strap
Nylon buckle

Composition :
SECURITY WRIST FASTENER
Combination of polyurethan foam and jersey/ cotton polyester fabric 100 % polyester velcro strap
100 % polyester strap

Composition :
ANKLE FASTENER
Combination of polyurethan foam and jersey/ cotton polyester fabric 100 % polyester fur
100 % polyester velcro strap
100 % polyester strap
Nylon buckle


PADDED WRIST FASTENER

LIGHT WRISTS SUPPORT FORR LOW-RESTLESSNESS PATIENTS

Set the wrist, on the fur side, in the bracelet.
With the help of the velcro, close the bracelet.
Wrap the strap around the object serving for fastening the fastener wrist.
Adjust the length of the belt with the nylon loop.

SECURITY WRIST FASTENER

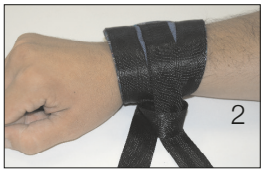FOR A STRONGER CONTENTION FOR VERY AGITATED PATIENTS.

Route the strap located at the sharp end of the bracelet into the slot (1).
Get the wrist of the person through this loop (2).
Fit to immobilize the wrist and close the bracelet using the Velcro, then move the strap through the loop (2).
Tie the straps around the object serving as fixing the security wrist fastener.

ANKLE FASTENER

MAINTAINS THE PATIENT'S ANKLES, WHILE ALLOWING FREEDOM OF MOVEMENT.

Set each ankle, on the fur side, in each bracelet.
Using the velcro, close each bracelet.
Route the strap around the object serving for fastening the ankle fastener.
Adjust the length of the belt with the nylon loop.

End-of-life management