

HAMMOCK SLING WITHOUT CROSSING

Avoid hip constriction.

The perforated hammock shape prevents the sling from being placed under the patient's buttocks, making donning and doffing more comfortable for the patient and easier for the carer.

The foam leg supports on this jersey and comfort mesh model reduce the pressure exerted on the thighs.

This sling prevents constriction of the hips and is suitable for short amputees. Each leg is held independently.

- With foam pads
- Available in comfort jersey and comfort mesh
- Sizes S, M, L, XL.


HOOK ADJUSTMENT



The multiple hooks system allows the patient's position to be varied.


For optimum comfort in the seated position, we recommend adjusting the sling as short as possible on the shoulders (green setting) and as long as possible on the legs (black setting).

SLING POSITIONNING

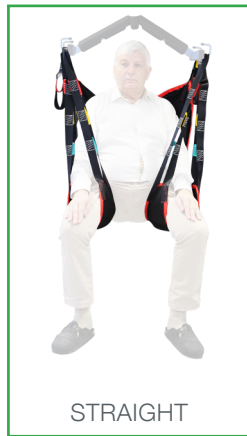

STRAIGHT

Recommended choice

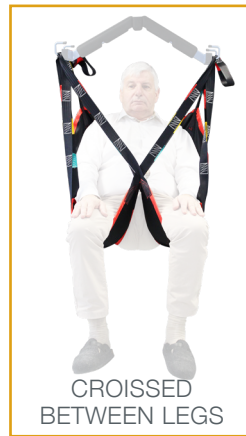

CROSSED BETWEEN LEGS

Second option

HEADREST (OPTION)

The headrest with clips (TAAT-PF) and removable headrest (TAMOV-PF) are available as options.



CONFORT STYLE BIEN-ÊTRE
FABRICATION FRANÇAISE



REF.	SUCPEHA2-PF (jersey comfort) SUCBAHA2-PF (mesh comfort)
OPTION	TAAT-PF (theadband with fastener) TAMOV-PF (removable headrest)