

PRODUCT SHEET

STAND-UP LIFT **WAY UP'GO COMPACT**

Descriptif :



- Improved handling and rolling thanks to a central wheel under the patient
- Motorized assistance for effortless movement under load
- Platform cushioning for improved threshold clearance
- Can be folded, cut and stored upright
- Steel frame
- Laser-cut and bent sheet metal, tubes and profiles
- Baked epoxy paint
- Ball-bearing wheels
- Total weight compact version: 32,5 kg
- Upper part weight : 17,40 kg (with one axle and two wheels)
- Weight of base plate : 15,10 kg (with other axle and two wheels)
- Maximum load : 180 kg
- Product life 8 years: 8 years* (excluding electrical part)
- * Subject to regular maintenance
- Lithium battery : 2-year warranty

- Medical devices (regulations (UE) 2017/745)
- Electromagnetic compatibility
- Compliance with the following standards :
 - NF EN 60601-1 : 2007 + A1 : 2014 + A2 : 2021
 - NF EN 60601-1-2 : 2016 + A1 : 2021
 - NF EN ISO 10535 : 2021
 - NF EN ISO 15223-1 : 2021
- 24V low-voltage DC motors
- Powerful 24 V/120 VA
- Maximum thrust force : 8 000 N
- Electronic overload protection
- Protection class : II
- Type of protection : IP X4
- Contactor limit switches
- ABS case
- Remote control with coiled cable, low voltage 24V
- Overload protection by thermo-monitors
- Noise level : 60 dB at 1 metre

Motor specifications



- Control box with removable high-capacity battery
- Graduated battery charge indication
- Autonomy indicator light
- Mains charger 220V-24V
- Emergency stop button
- 4-function remote control IPX4 protection (up/down, forward/backward))
- Emergency lowering : electric and manual B D
- Mains charger connection for battery charging

Dimensional characteristics

- A.** Total length : 91 cm
- B.** Minimum useful height : 87 cm
- C.** Max. useful height : 124.5 cm
- D.** Total height (excluding remote control) : 91 cm
- E.** Frame height : 8 cm
- F.** Minimum base width : 39 cm
- G.** Max. baseplate width: 58 cm
- H.** Turning diameter : 100 cm
- I.** Ground clearance : 4 cm

In upright folded storage (see photo opposite) :

- A.** Width: 59 cm
- B.** Depth : 49 cm
- C.** Height : 110 cm

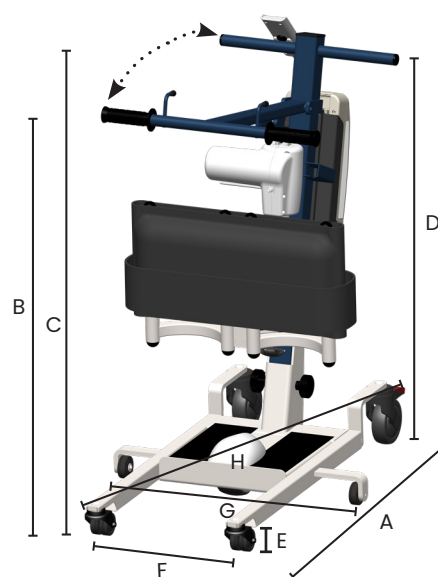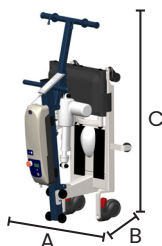

Non-contractual visual

Packaging : 1000x600x600 mm (stackable up to C three products)

- Pallet 800x1200 mm : up to 3 products
- Pallet 1000x1200 mm : up to 6 products